



DynaFlo 8000 Dynamic Pressure Relief Mattress Replacement System



This 8" replacement mattress with cell-in-cell designed air bladders offers a low air loss therapy. The patented cell-in-cell design allows the inflation and deflation of air cells to be achieved in the shortest amount of time while providing the good pressure dispersion. The inner cells also offer a sufficient weight support during a power failure. The microprocessor pressure sensor technology provides smart pressure auto-checking feature to maximize patients comfort. This system is suitable for treatment of pressure sore stage I, II & III prevention with maximum weight capacity of 160kg.

Key Features:

- ◆ Microprocessor pressure sensing technology to produce low pressure support.
- ◆ Visible and audible alarms with Alarm Mute feature.
- ◆ Transport mode by interconnecting quick connectors.
- ◆ External, large air filter intake for increased life expectancy and easy maintenance.
- ◆ Static function can suspend the alternating mode.
- ◆ Spring loaded tube to prevent any kinks for smoother airflow.
- ◆ Elegant zipper-all-around cover for easy cleaning and maintenance.
- ◆ Patented Cell-in-Cell bladder to provide adequate pressure support during power failure.
- ◆ Low air loss feature helps to prevent heat and moisture build-up.
- ◆ Vapour permeable, water resistant cover with non-woven quilted material underneath.
- ◆ Slip-resistant mat and securing straps on the mattress base for added safety.



Microprocessor



Static Button



Large Air Filter Intake



Quick Connectors




Air Bladders with Low Air Loss



Zipper-All-Around Cover



Patented Cell-in-Cell bladder's design

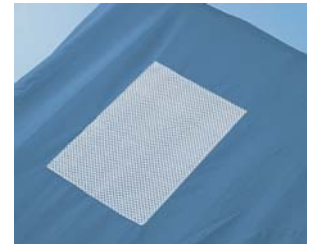

Slip-Resistant Mat

Specifications

DynaFLO 8000

PUMP	Electrical	AC120V/60Hz, <1A AC230V/50Hz, <1A	
	Air Output (lpm)	> 8	
	Pressure Range (mmHg)	10 - 40	
	Cycle Time (min)	10 min	
	Dimension (L x W x H)	30 x 10.5 x 22.5 cm / 11.8" x 4.1" x 8.9"	
	Weight	3.4 kg / 7.5 lb	
MATTRESS	Arrangement	20 x 8" cells with cell-in-cell design	
	Material	Cover	Nylon / PU , dark blue
		Cells	Nylon / PU , light blue
		Base	Nylon / PU , dark gray
	Dimension (L x W x H)	200 x 90 x 20.3 cm or 78.7" x 35.4" x 8"	
Weight	10 kg / 22 lb		

All specifications are subject to change without notice.



Private label is welcome



APEX MEDICAL CORP.

9, Min Sheng St. Tu-Cheng,
 Taipei County, 236, Taiwan, R.O.C.
 TEL: 886-2-22683100
 FAX: 886-2-22686525
<http://www.apexmedical.com>