

# NeoBLOW 8200

## True Low Air Loss & Dynamic Flotation Mattress Replacement System

### NeoBLOW 8200

The NeoBLOW 8200 True Low Air Loss Replacement System uses a high volume airflow blower and a numerous air loss holes on bladders for achieving the true low air loss therapy to provide excellent ventilation and keep patients dry from humidity. The 10 inches high mattress, which consists of 8 inches high air bladders at the top and a layer of 2" foam at the bottom, provides maximum comfort for patients and prevents from bottoming out. The pulsation function on the control unit provides patients the massage action to stimulate the blood circulation and reduce the edema. The system supports patients up to 180 kg and is suitable for treatment of pressure sore stage IV.


### Key Features:

- ◆ True low air loss and dynamic flotation therapy.
- ◆ Smart sensor technology to produce the lowest pressure support.
- ◆ AutoFirm function to inflate mattress quickly and facilitate patient to access.
- ◆ Panel Lock function prevents persons other than caretakers from changing the setting.
- ◆ Pulsation feature provides gentle massage to stimulate the blood circulation.
- ◆ 3 different pressure zones: head, body & leg with 6 comfort level adjustments.
- ◆ Easy operating with easy recognized graphic icons.
- ◆ Vapour permeable, water resistant cover.
- ◆ 2" foam base provide the maximum safety to prevent bottoming out during power failure.
- ◆ Slip-resistant mat and securing straps on the mattress base for added safety.




#### NeoBLOW Pump

Easy operating with easy recognized graphic icons.

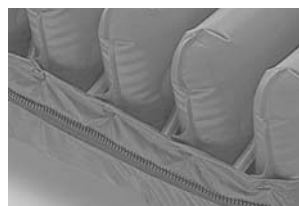

#### Zipper-All-Around Cover

Elegant looking finish with mat tress set. Cover can be taken out fir cleaning or replacing.

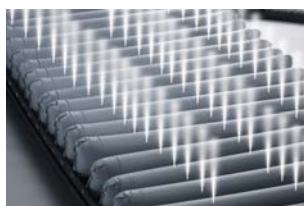

#### True Low Air Loss

A high volume airflow blower pro vide excellent ventilation and keep patients dry from humidity.



#### Large Air Filter Intake

Easy to change, washable filter.

### Specifications

#### NeoBLOW 8200

		NeoBLOW 8200	
PUMP	Electrical	AC120V/60Hz, <5A AC230V/50Hz, <5A	
	Air Output (lpm)	Blower Technology	
	Pressure Range (mmHg)	Adjustable 6 levels	
	Cycle Time (min)	N/A	
	Dimension (L x W x H)	37x27x13cm	
	Weight	4.5 kg	
MATTRESS	Arrangement	20x8" + 2" foam pocket	
	Material	Cover	Nylon/PU
		Cells	Nylon/PU
		Base	Nylon/PU
	Dimension (L x W x H)	200x90x20.3cm	
Weight	7.0Kg		



All specifications are subject to change without notice.

Private label is welcome



## APEX MEDICAL CORP.

9, Min Sheng St. Tu-Cheng,  
Taipei County, 236, Taiwan, R.O.C.  
TEL: 886-2-22683100  
FAX: 886-2-22686525  
<http://www.apexmedical.com>